

EVIDENCE MATTERS

RESEARCH BULLETIN

Stryker® MicroFX™ OCD (OsteoChondral Drill) Results in Less Bone Compaction Around Microfracture Holes¹



TOP-LEVEL SUMMARY

The density of the bone was compared around microfracture holes created with the Stryker MicroFX OCD and the Linvatec awl. **The Stryker MicroFX OCD drill resulted in significantly less bone compaction immediately around the microfracture hole.**¹

METHODS


Microfracture holes were created in the condyles of six goats using either a 45 degree MicroFX OCD or a 45 degree awl. After 3 days the knees were harvested and the density of the bone in a cylindrical area from 0 to 0.5 mm adjacent to the microfracture hole was evaluated using microCT (uCT 80; SCANCO Medical, Bassersdorf, Switzerland). Significant differences were assumed at $p \leq 0.05$.

RESULTS

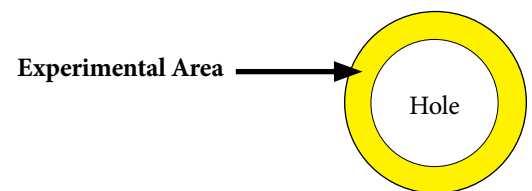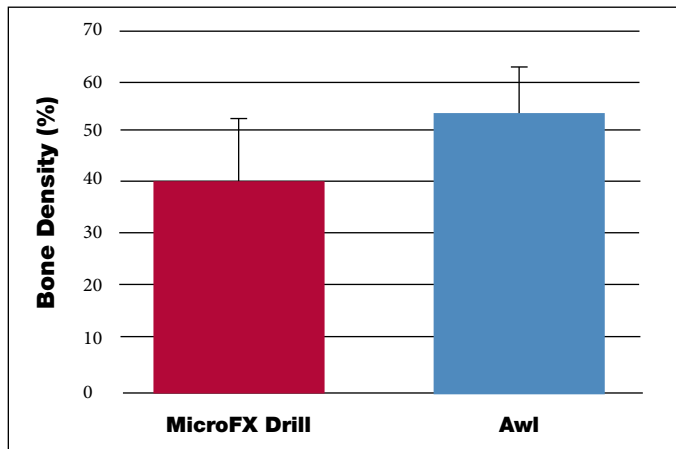


The bone around microfracture holes created with the drill was significantly less dense and had a consistent density when compared to the bone around holes created with the awl.

The bone around the microfracture holes created by the awl was very dense on one side - the side away from the direction of impact (shown as a yellow arrow) and was less dense on the opposite side.

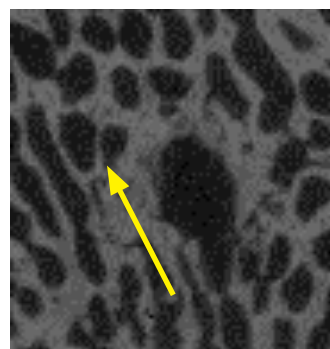

Awl

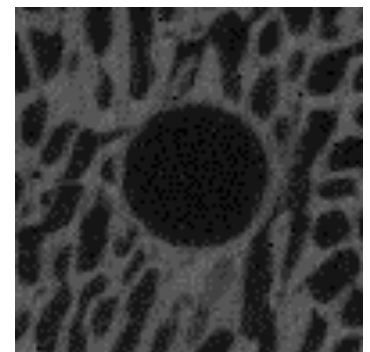

MicroFX OCD Drill

CLINICAL RELEVANCE

The Stryker MicroFX OCD drill resulted in significantly less bone compaction immediately around the microfracture hole providing more open channels for blood flow into the defect. Since the objective of creating microfracture holes is to encourage the flow of blood (with healing elements) into the defect in order to facilitate healing of chondral defects,² use of the OCD drill may result in excellent blood flow.

Joint Replacements

Trauma, Extremities & Deformities

Craniofacial

Spine

Biologics

Surgical Products

Neuro & ENT

Interventional Spine

Navigation

Endoscopy

Communications

Imaging

Patient Care & Handling Equipment

EMS Equipment

325 Corporate Drive
Mahwah, NJ 07430
t: 201.831.5000

www.stryker.com

References

1. Technical Report #RD13-019 Rev. 1
2. Steadman JR, Rodkey WG, Briggs KK "Microfracture:
Its History and Experience of the Developing Surgeon" Cartilage 1(2) 78-86, 2010.

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: MicroFX, Stryker. All other trademarks are trademarks of their respective owners or holders.